

EXHIBIT 2



PROFESSIONAL HEALTH

CONSUMER HEALTH

ENGLISH SEARCH

TECHNOLOGY

PRODUCT

SOLUTION

CLINICAL EVIDENCE

SUPPORT

Home / Root Platform / Next Generation SedLine® Brain Function Monitoring



Next Generation SedLine® Brain Function Monitoring

A More Complete Picture Starts with More Complete Data

Root® with Next Generation SedLine Brain Function Monitoring helps clinicians monitor the state of the brain under anesthesia with bilateral data acquisition and processing of electroencephalogram (EEG) signals.

Next Generation SedLine Features:


- › Four simultaneous channels of frontal EEG waveforms
- › An enhanced Patient State Index (PSi), a processed EEG parameter related to the effect of anesthetic agents
- › A Density Spectral Array (DSA) display, which contains left and right spectrograms representing the power of the EEG on both sides of the brain
- › An optional Multitaper DSA, which may enhance visibility of EEG features

Improved Patient State Index (PSi)

Next Generation SedLine features an enhanced signal processing engine which provides an enhanced Patient State Index (PSi), a processed EEG parameter related to the effect of anesthetic agents.

EEG experts scored the improvement in PSi performance between the original SedLine PSi and Next Generation SedLine PSi. A 25% average improvement in Next Generation PSi performance was found.

Expert Scoring of Next Generation SedLine¹

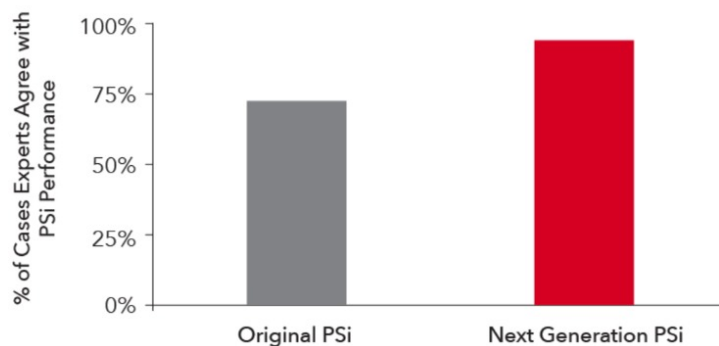

To evaluate the performance of Original PSi and Next Generation PSi, independent EEG experts reviewed cases with both Original PSi and Next Generation PSi (blinded to the version), along with additional clinical information (MOAAS scores, EEG waveforms, drug doses, and vital signs). Compared to the expert-assessed anesthetic depth, an error was defined as a case when expert assessment of PSi was 'Low' or 'High' and success was defined as a case when the expert assessment of PSi was 'Good'.


Parallel Signal Processing Engines

EMG is a common confounding factor that can interfere with EEG signals used in brain function monitoring.² Next Generation SedLine utilizes Masimo's Parallel Signal

Adaptive Signal Processing with Band-Independent Features


Power across all frequency bands decreases with age.³ When computing PSi, Next Generation SedLine uses adaptive signal processing with band-independent features

Processing Engines to extract a clearer EEG signal for computing PSI.

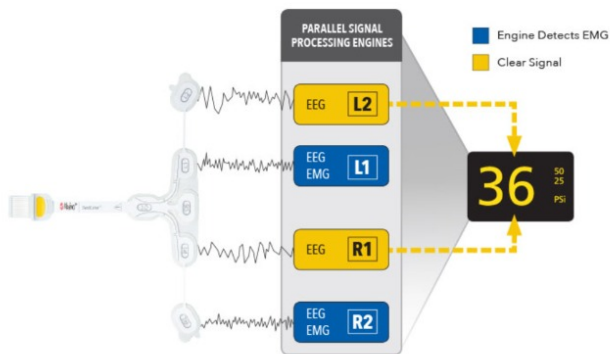


This image captures a moment when Next Generation SedLine detects EMG in the two engines depicted.

to search for EEG features across many frequency bands.

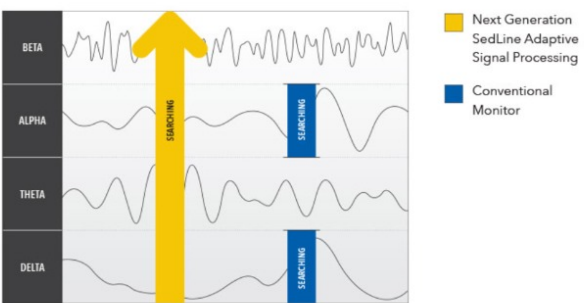

Multitaper Density Spectral Array (DSA)

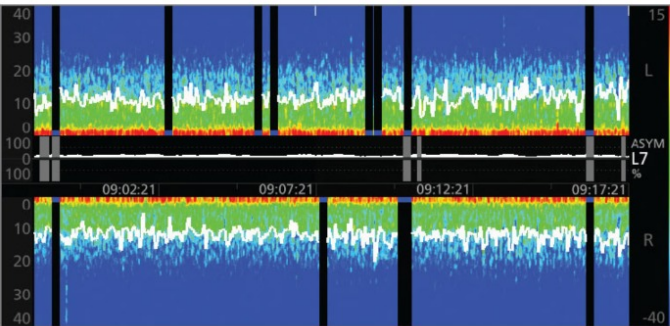

Next Generation SedLine offers clinicians the flexibility of choosing to display either an enhanced Multitaper Density Spectral Array (DSA) or a standard Hanning DSA. The DSA contains left and right spectrograms representing the power of the EEG on both sides of the brain.

When using a Multitaper DSA, EEG data are transformed into the frequency domain, which may provide a better display of EEG features.

A More Complete Picture of the Brain

The Next Generation SedLine module easily plugs into the Root patient monitoring platform via Masimo Open Connect® (MOC-9®) ports. Root's customizable, easily-interpretable display offers multiple views of brain monitoring information expanding visibility in the operating room and intensive care unit.



Next Generation SedLine can be used simultaneously with O3® Regional Oximetry on the Root platform for a more complete picture of the brain.




RD SedLine™ EEG Sensor




- › Four active EEG leads collect data in the frontal lobe
- › Soft foam pads for comfortable application on a patient's forehead
- › Adhesive anchors ensure secure sensor placement for optimal signal quality
- › Pre-filled gel electrodes help streamline sensor application workflows
- › Application graphics for O3 regional oximetry sensor placement simplify simultaneous application of both monitoring technologies

Brain Monitoring Solutions


Explore the Root Platform and Brain Monitoring Sensors & Cables

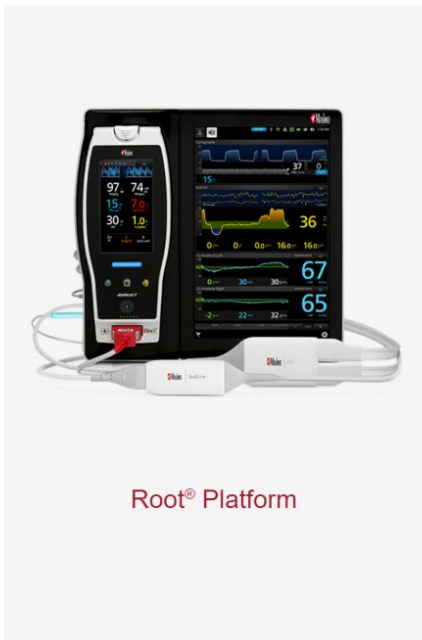

Root® Platform

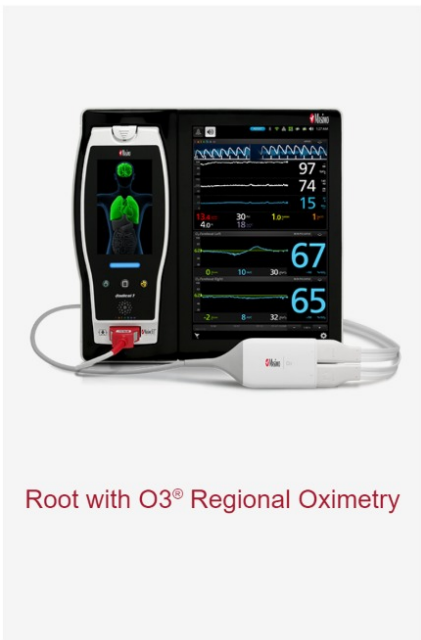

Root with O3® Regional Oximetry

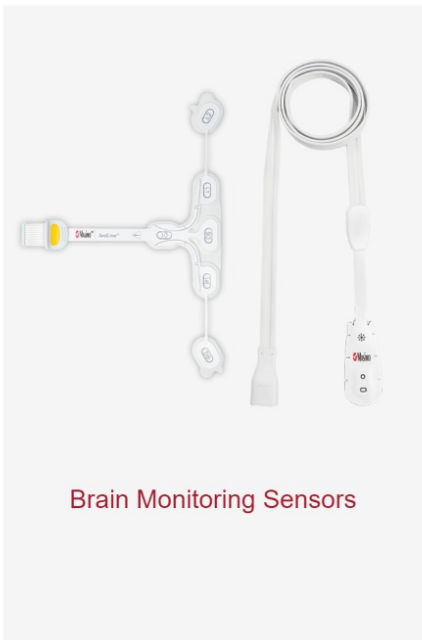















Brain Monitoring Sensors

REQUEST MORE INFORMATION

References:

¹ Analyse rétrospective des données cliniques.
² Lobo, Francisco A., and Stefan Schraag. Limitations of anaesthesia depth monitoring. *Current Opinion in Anesthesiology*. 24, no. 6 (2011): 657-664.
³ Purdon P et al. *Brit J of Anaesth*. 10.1093 46-57.

RESOURCES

	<div>Next Generation SedLine</div> <div>Brochure</div>	
	<div>Root with SedLine Brain Function Monitoring</div> <div>Product Information</div>	
	<div>Perioperative Care Solutions</div> <div>Brochure</div>	
	<div>RD SedLine Sensor</div> <div>Product Information</div>	
	<div>Guide SedLine</div> <div>Quick Reference</div>	
	<div>Next Generation SedLine</div> <div>Video</div>	
	<div>Brain Monitoring - DSA Interpretation</div>	

RD SedLine sensor is not licensed for sale in Canada.

Caution: Federal (USA) law restricts this device to sale by or on the order of a physician. See instructions for use for full prescribing information, including indications, contraindications, warnings, and precautions.





Masimo

+ Improve Life™

+ Improve patient outcomes and reduce cost of care®

+ Take noninvasive monitoring to new sites and applications™

Follow Us



ABOUT US

Contact Us

Guiding Principles

Corporate Brochure

Sustainability Report

Company Evolution

Investor Relations

Press Room

Media Room

Careers

Newsletter sign-up

FOR HEALTHCARE

Technology

Product


Solutions


Clinical Evidence


OEM Solutions


Sensors and Cannulas


Support


Resources 

FOR CONSUMERS 

Products 

Experience 

Education 

Support 

LEGAL

Compliance

Codes of Conduct

Patents

Terms

PRIVACY

Privacy Notice

CA Privacy

Cookie Notice

Cookie Preferences

SITEMAP